

Physiotherapy Equipments

▶ Heat & Cold Therapy Equipments

We have quality products in Heat & Cold Therapy equipments. All our models are efficient and display quality in Physiotherapy Equipments.



IMI - 2515: PARAFFIN WAX BATH (Arm, Hand & Foot)



Paraffin treatments are used for the Symptomatic relief of pain and stiffness due to Arthritis, Bursitis, Tendonitis, Muscle strains or Sports Injuries. The low temperature paraffin permits the immersion of the limb directly into the wax-bath or to apply with brush.

Inner size of S.S. tank:
50cm x 35cm x 25cm.

[\[Request Quotation / More info...\]](#) 🏠

SPECIFICATIONS:

*Construction: Double walled effectively insulated. Wax tank of 20 gauge Stainless steel sheet. Fitted with laminated wooden rim & anodized Aluminium cover Mounted on 5cm dia. Four casters.

*Heaters: Two, 1kw each, special wax Heaters.

*Heat Control: Thermostat to control temperture from 30 to 110 deg.Cen.

*Finish: Externally finished in Oven Baked.

*Accessories: 6 Kg. Paraffin Wax, Power Plug, One brush for wax application & aluminium sheet mesh for heaters.

IMI - 2516: PARAFFIN WAX BATH (Arm & Hand)



Internal wax tank size : 40cm.x 20cm.x 20cm. deep, is made out of stainless steel sheet.The top has a metallic cover and the edges of the bath has a laminated covered wooden rim, all around the top. Wax bath moves on 50mm. dia casters.Works on 220V A.C. Provided with one special wax immersion heater and thermostate to control temperature from 30 to 110 deg.Cent.

Finished in oven baked enamel. Comes complete with three core wire,power pin plug & 3kg wax.

[\[Request Quotation / More info...\]](#)

IMI - 2518: PHYSIO-WAX

No picture available

Special quality formulated wax with low melting point i.e. about 52 deg.'C' for use in paraffin wax bath.

[\[Request Quotation / More info...\]](#)

IMI - 2521: MOIST HEAT THERAPY UNIT (Table Model)



Made out of heavy gauge stainless steel sheet double walled well insulated in between.

Overall size 38cmx 30cmx 45cm high.

Fitted with 1000 watts immersion heater, pilot light and thermostat for heat control. Comes with 4 standard size 25cmx30cm steam packs containing special heat retaining silica gel crystals. Polished finish inside and special paint finish outside.

[\[Request Quotation / More info...\]](#)

IMI - 2522: MOIST HEAT THERAPY UNIT (Junior, Mobile)



Made out of heavy gauge stainless steel sheet double walled, fully insulated in between top and sides.

Overall size 45cm x30cm x75cm high.

Fitted with 2000 watts immersion heater, pilot lights and thermostat for heat control. A drain valve is provided for emptying. Mounted on 75mm dia rubber castors for easy mobility.

Comes with 8 Steam packs (five of standard size 25cm x30cm, one of half size 12cm X 30cm, one of contour shape cervical pack size 16cm X 60cm and one full back size 25cmx50cm, containing special heat retaining silica gel crystals.

Polished finish inside & spc. paint finish outside.

[\[Request Quotation / More info...\]](#) 🏠

IMI - 2523: MOIST HEAT THERAPY UNIT (Senior, Mobile)



The specifically designed heating unit automatically maintain Hot Packs in water at proper temperature with no plumbing necessary.

[\[Request Quotation / More info...\]](#) 🏠

SPECIFICATIONS:

*Construction: Unit's outer body made of 18 gauge & inner chamber of 20 gauge S.S.Sheet double walled,well insulated, mounted on 75mm dia. 4 casters & drain valve.

*Dimensions: Over all, approx.60cm x40cm x80cm.

*Heat Control: Thermostate to control Temperature from 30 to 110 deg.Cen.with 2 Pilot Lamps.

*Heater: 2 kw. specifically designed Heater.


*Hot Packs: Total 12 of different sizes as

- a) Six, Standard size (25x30cm)
- b) Two, Half Size (12x30cm)
- c) Two, Neck Contour (60x15cm)
- d) Two, Full Back (25x50cm)

*Accessories: Two forceps to get-out Hot-packs from hot water & Power Pin Plug.

*Finish: Polished finish inside.

:: STEAM / HOT PACKS

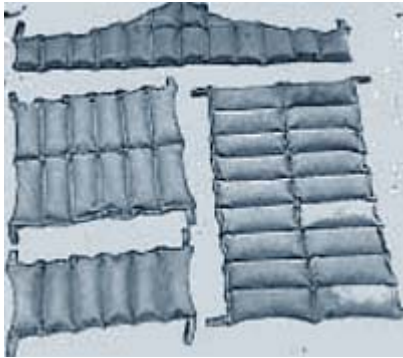

Steam / Hot packs are made of strong water absorbing fabric containing special silica gel crystals which absorb many times its own volume of water & retain for pretty long time. Packs can be applied again & again for moist heat.

Standard size 25cm x 30cm.

[\[Request Quotation / More info...\]](#) 🏠

IMI - 2525: STEAM / HOT PACKS (Standard Size: 25 x 30 cm.)

IMI - 2526: STEAM / HOT PACKS (Half Size: 12 x 30 cm.)

IMI - 2527: STEAM / HOT PACKS (Contour Shape: 16 x 60 cm.)

IMI - 2528: STEAM / HOT PACKS (Full Back Size: 25 x 50 cm.)

IMI - 2536: IMICO COLD THERAPY UNIT (with Fan)



Made out of heavy gauge stainless steel sheet. Tank is fully insulated in between top and sides. Fitted with heavy duty cooling unit & thermostate for temperature control. A drain valve is provided for emptying. Mounted on 75mm dia. ball bearing castors for easy mobility. Comes with SIX standard size Cold packs containing special gel.

Finished bright inside and special paint finish outside. It completely simplifies the problem of Cold Therapy applications in physcial-therapy departments.

[\[Request Quotation / More info...\]](#) 

IMI - 2538: CONTRAST BATH (Hot & Cold therapy unit)



Suitable for treatment of Arms & Legs. Two Stainles Steel heavy gauge tanks with inside measurements of 50cm x35cm x 38cm deep each are mounted on 10cm dia heavey ball bearing castors. One tank is fitted with special heating elements, the other with heavy duty cooling compressure unit. Each tank provided with Digital Thermometer, Thermostat for Temperatur Control, Pilot Lights & Water Outlet.

Both tanks are well insulated and placed side by side in a heavy tubular steel frame covered with painted steel sheet.

Works on 220 V. A.C. UNIT COMES WITH 5 kv. VOLTAGE STABILISER.

[\[Request Quotation / More info...\]](#) 